





Constant climate chamber HPP260

Our constant climate chamber is tailored to environmental simulation, material testing and stability testing in accordance with the ICH guidelines.



On this page, you can find all the essential technical data on the Memmert stability chamber HPP. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at myAtmoSAFE@memmert.com.

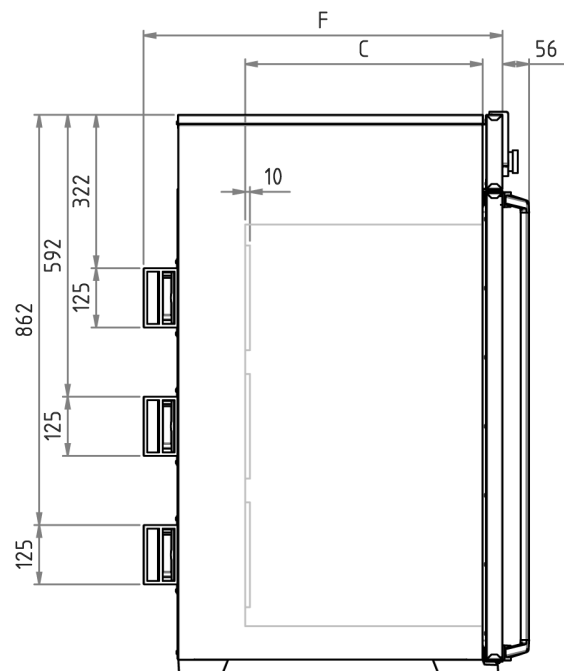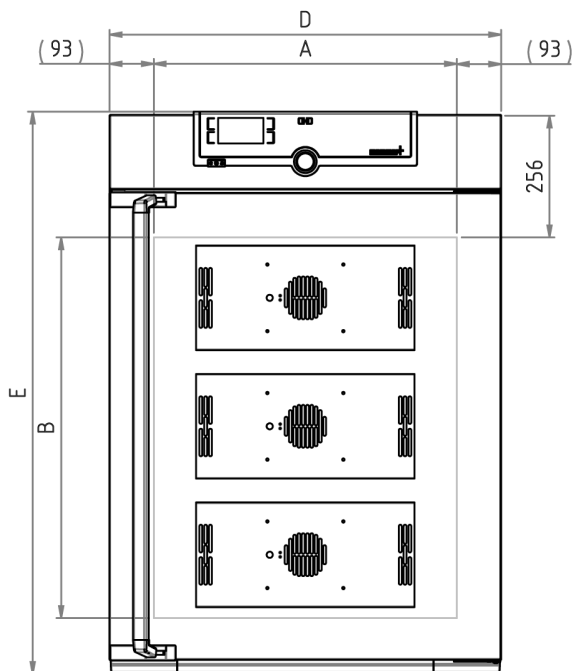

Temperature

Working temperature range	without light: from 0°C to +70°C
Working temperature range	with light: from +15°C to +40°C
Setting accuracy temperature	0.1°C
Temperature	2 Pt100 sensors Class A in 4-wire-circuit, mutually monitoring and taking over the performance at the same temperature value

Humidity

Setting range humidity	10 - 90 % rh
Setting accuracy humidity	0.5 %
Humidity	humidity supply with distilled water from external tank by self-priming pump
Humidification	humidification by hot steam generator
Dehumidification	dehumidification by cold trap using the Peltier technology

Control technology

ControlCOCKPIT	TwinDISPLAY. Adaptive multifunctional digital PID-microprocessor controller with 2 high-definition TFT-colour displays.
Language setting	German, English, Spanish, French, Polish, Czech, Hungarian
adjustable parameters	temperature (Celsius or Fahrenheit), relative humidity, programme time, time zones, summertime/wintertime
Timer	Digital backwards counter with target time setting, adjustable from 1 minute to 99 days
Function HeatBALANCE	adapting the distribution of the heating performance of the upper and lower heating circuit from -50 % to +50 %
Function SetpointWAIT	the process time does not start until the set temperature is reached
Calibration	three freely selectable temperature values, 2-point calibration for humidity: 20 and 90 % rh

Communication

Interface	Ethernet LAN, USB
Documentation	programme stored in case of power failure
Programming	AtmoCONTROL software on a USB stick for programming, managing and transferring programmes via Ethernet interface or USB port

Safety

Temperature control	over- and undertemperature monitor TWW, protection class 3.3 or adjustable temperature limiter TWB, protection class 2, selectable on display
AutoSAFETY	additionally integrated over- and undertemperature protection "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating function is switched off in case of overtemperature, cooling function in case of undertemperature
Autodiagnostic system	integral fault diagnostics for temperature and humidity control
Alarm	visual and acoustic

Heating concept

Peltier	energy-saving Peltier heating-/cooling system integrated in the rear (heat pump principle)
----------------	--

Standard equipment

Standard accessories	Water tank including connection hose
-----------------------------	--------------------------------------

Works calibration certificate	standard values +10°C u. +37°C, 60 % rh bei +30°C
--------------------------------------	---

Door	fully insulated stainless steel door with 2-point locking (compression door lock)
-------------	---

Door	inner glass door
-------------	------------------

Internals	2 stainless steel grids
------------------	-------------------------

Stainless steel interior

Volume	256 l
---------------	-------

Dimensions W x H x D in mm	$w_{(A)} \times h_{(B)} \times d_{(C)}$: 640 x 800 x 500 mm
-----------------------------------	--

Max. number of grids/shelves	9
-------------------------------------	---

Max. loading of chamber	200 kg
--------------------------------	--------

Max. loading per grid/shelf	20 kg
------------------------------------	-------

Textured stainless steel casing

Dimensions	$w_{(D)} \times h_{(E)} \times d_{(F)}$: 824 x 1183 x 774 mm
-------------------	---

Housing	rear zinc-plated steel
----------------	------------------------

Electrical data

Voltage	230 V, 50/60 Hz
Electrical load	approx. 920 W

Voltage	115 V, 50/60 Hz
Electrical load	approx. 920 W

Packing/shipping data

the appliances must be transported upright

Customs tariff number	8419 8998
------------------------------	-----------

Country of origin	Federal Republic of Germany
--------------------------	-----------------------------

WEEE-Reg.-No.	DE 66812464
----------------------	-------------


Dimensions approx incl. carton	B x H x T: 930 x 1380 x 930 mm
---------------------------------------	--------------------------------

Net weight	approx. 122 kg
-------------------	----------------

Gross weight carton	approx. 173 kg
----------------------------	----------------

Temperature-humidity working range HPP

Not all climate chambers are the same. The humidity content of the chamber load, the ambient conditions and the respective temperature-humidity working range are decisive factors in the selection of the right appliance. In the adjacent diagram, you can see the possible temperature/humidity combinations for our constant climate chamber HPP.

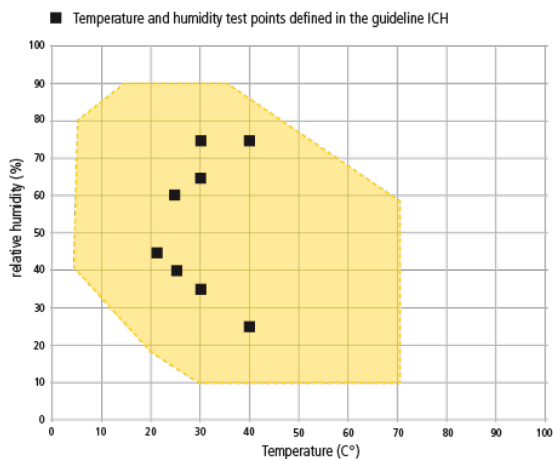

Temperature-humidity working range HPP

Standard units are safety-approved and bear the test marks

