




Environmental test chambers

TTC256

Memmert environmental test chambers guarantee the perfect atmosphere for climate and temperature tests



On this page, you can find all the essential technical data on the Memmert environmental chamber CTC/TTC for climate testing and temperature testing. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at myAtmoSAFE@memmert.com.

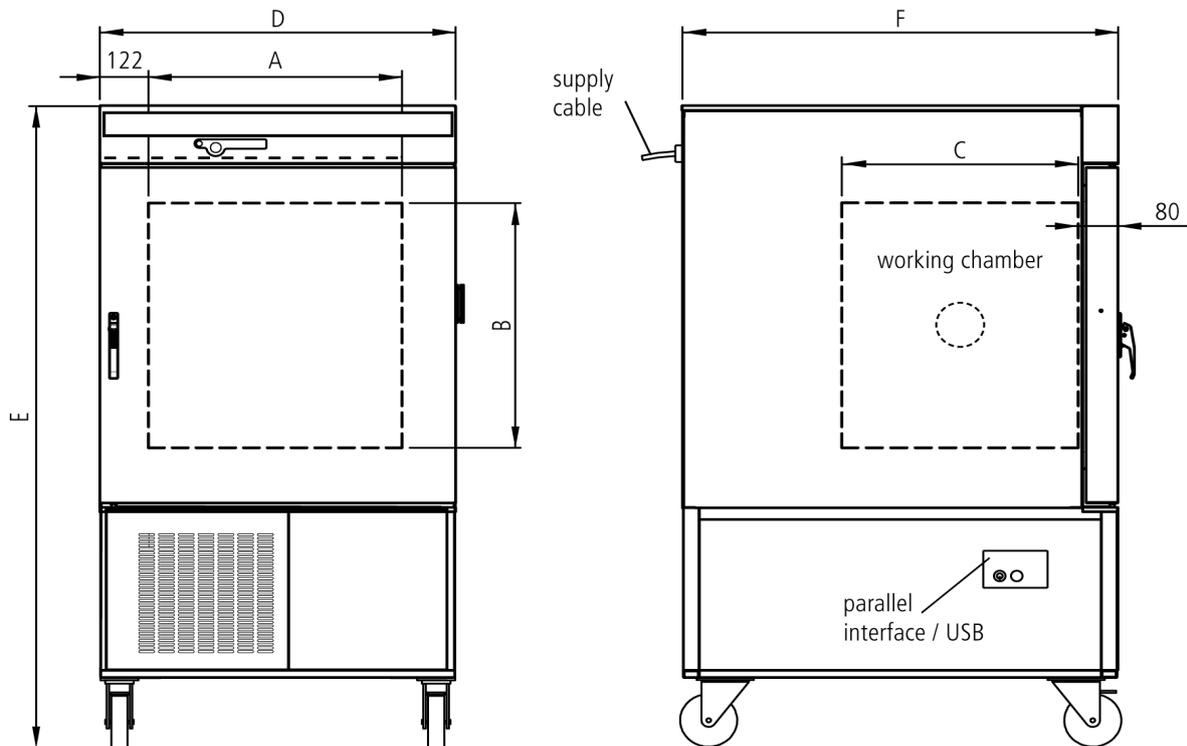

Temperature

Temperature range	without humidity from -42°C up to +190°C
Temperature distribution (spatial)	+/- 0.5 up to 2K
Temperature	2 Pt100 sensors Class A in 4-wire-circuit, mutually monitoring and taking over the performance at the same temperature value
Display	resolution of display for setpoint values 0.1°C up to 99.9°C, 0.5°C from 100°C and for actual values 0.1°C (LED)

Humidity

Dehumidification	2 x 10 l tanks as condensate collector
-------------------------	--

Control technology

Controller	Electronic microprocessor temperature controller with auto-diagnostic system
Controller	digital display of all set parameters, such as temperature, weekdays, time, fan speed, air flap position, programme status and set-up values - language to be chosen via set-up
Timer	integrated timer for tempering profiles of up to 40 ramps each, each segment adjustable from 1 min. to 999 hrs.
Calibration	three freely selectable temperature values

Ventilation

Fan	high-performance air fan in working chamber
Fan control	depending on operation status automatic adaption of fan speed resp. manual adjustment from 10 to 100 %

Communication

Interface USB	USB-interface incl. Memmert software "Celsius" for programming and documentation
Interface Printer	parallel printer interface (incl. real time clock with date function) for all PCL3-compatible ink jet printers for GLP-conforming documentation
Documentation	integrated ring memory as data logger for GLP-conforming long-term documentation of all relevant parameters - 1024 kB
Documentation	programme stored in case of power failure
Programming	chip-card control incl. 1 MEMoryCard XL with 32 kB storage capacity (max. 40 ramps)

Safety

Temperature control	mechanical temperature limiter TB, protection class 1 according to DIN 12880 to switch off the heating approx. 10°C above nominal temperature
Temperature control	over- and undertemperature monitor TWW, protection class 3.3 or adjustable temperature limiter TWB, protection class 2, adjustable at controller
AutoSAFETY	additionally integrated over- and undertemperature protection "ASF", automatically following the setpoint value at a preset tolerance range, alarm in case of over- or undertemperature, heating is switched off in case of overtemperature, compressor in case of undertemperature
Autodiagnostic system	integral fault diagnostics for temperature control
Alarm	with visual and acoustic alarm in case of over-/under temperature and open door

Heating concept

Ring heater high-performance ring heaters with optimised air circulation

Door door heating to avoid condensation

Cooling twin compressor; chlorine-free refrigerant R404A

Standard equipment

Works calibration certificate standard values +30°C u. 60 % rh

Door fully insulated stainless steel door with double-locking and 4-point adjustment, heated

Interior easy-to-clean stainless steel interior, materil 1.4301 (ASTM 304), hermetically welded

Internals 1 stainless steel grid

Stainless steel interior

Volume 256 l

Dimensions W x H x D in mm $w_{(A)} \times h_{(B)} \times d_{(C)}$: 640 x 670 x 597 mm

Max. number of grids/shelves 6

Max. loading per grid/shelf 25 kg

Textured stainless steel casing

Installation on lockable castors

Dimensions $w_{(D)} \times h_{(E)} \times d_{(F)}$: 898 x 1730 x 1100 mm

Housing rear zinc-plated steel

Electrical data

Electrical load 400 V, 3 phases, 50/60 Hz / approx. 7000 W

Packing/shipping data

the appliances must be transported upright

Customs tariff number 8419 8998

Country of origin Federal Republic of Germany

WEEE-Reg.-No. DE 66812464


Dimensions approx incl. carton B x H x T: 1020 x 1910 x 1310 mm

Net weight approx. 337 kg

Gross weight carton approx. 463 kg

Temperature-humidity working range CTC

Not all climate chambers are the same. The humidity content of the chamber load, the ambient conditions and the respective temperature-humidity working range are decisive factors in the selection of the right appliance. In the adjacent diagram, you can see the possible temperature/humidity combinations for our climatic test chamber CTC.

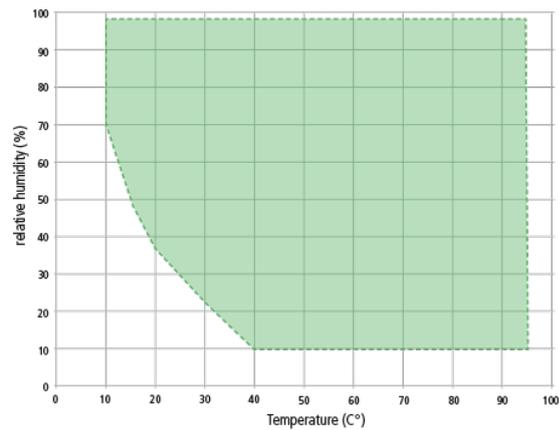

Temperature-humidity working range CTC

Standard units are safety-approved and bear the test marks

