




Waterbath

WNE 29

The combination of corrosion-resistant stainless steel, precise electronics and multiple temperature protection guarantees the highest level of safety in the laboratory.



On this page, you can find all the essential technical data on the Memmert water bath. Our customer relations team will be pleased to help if you want further information. If you should require a customised special solution, please contact our technical specialists at myAtmoSAFE@memmert.com.

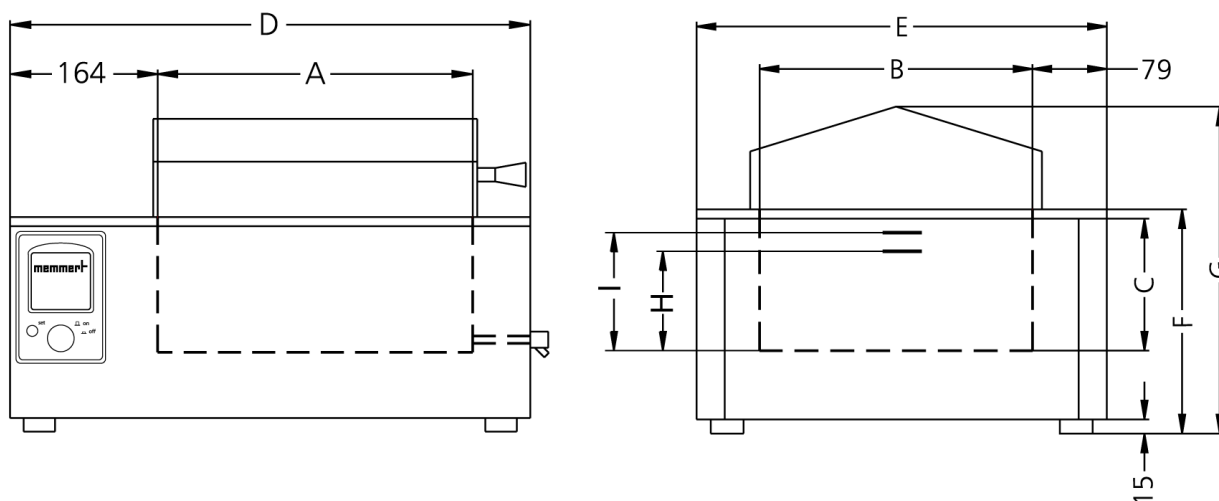

Temperature

Temperature range in °C	min. 5°C above ambient up to +95°C with additional boiling mode (+100°C)
resolution of display and setting accuracy	+0,1°C below 99,9°C, 1°C above 100°C

Control of standard components

Controller	digital display (LED) of all set parameters, such as temperature and alarm values (0,1°C resolution) as well as time values
Timer	digital timer from 1 min. up to 999 hours for: ON (continuous operation), DELAYED ON, HOLD or HOLD set-temperature dependent with guaranteed holding time
Calibration	on controller

Safety

Temperature control	2 Pt100 sensors Class A in 4-wire-circuit, mutually monitoring and taking over the performance at the same temperature value
Temperature control	mechanical temperature limiter TB protection class 1 switching the heating off at approx. 30°C above max. temperature of the bath
Temperature control	independently working, electronic overtemperature controller TWW protection class 3.1 or overtemperature limiter TWB protection class 2, adjustable in set-up menu by customer
Temperature control	in case of overtemperature due to failure, the heating is switched off at approx. 10°C above the set temperature (fixed value)
Autodiagnostic system	fuzzy-supported PID microprocessor controller with integrated autodiagnostic system with fault indication
Alarm	visual and acoustic alarm at programme end and as input acknowledgement as well as in case of low liquid level, heating is switched off automatically

Heating concept

Heating Baths	corrosion-proof large-area heating on three sides
----------------------	---

Stainless steel interior

Interior	easy-to-clean interior, made of stainless steel, reinforced by deep drawn ribbing, material no. 1.4301 (ASTM 304), laser-welded
Volume	29 l
Dimensions W x H x D in mm	$w_{(A)} \times h_{(B)} \times d_{(C)}$: 590 x 350 x 140 mm
Liquid level min.	(H) 105 mm
Liquid level max.	(I) 120 mm

Textured stainless steel casing

Dimensions	$w_{(D)} \times h_{(E)} \times d_{(F)}$: 818 x 516 x 238 mm
-------------------	--

Electrical data

Voltage	230 V, 50/60 Hz
----------------	-----------------

Electrical load	approx. 2400 W (during heating)
------------------------	---------------------------------

Packing/shipping data

the appliances must be transported upright

Customs tariff number	8419 8998
------------------------------	-----------

Country of origin	Federal Republic of Germany
--------------------------	-----------------------------

WEEE-Reg.-No.	DE 66812464
----------------------	-------------

Dimensions approx incl. carton	B x H x T: 910 x 610 x 400 mm
---------------------------------------	-------------------------------

Net weight	approx. 24 kg
-------------------	---------------

Gross weight carton	approx. 31 kg
----------------------------	---------------

Standard units are safety-approved and bear the test marks

